


INVENTED  
FOR  
CHILDREN



DELFI PEDIATRIC TOURNIQUETS

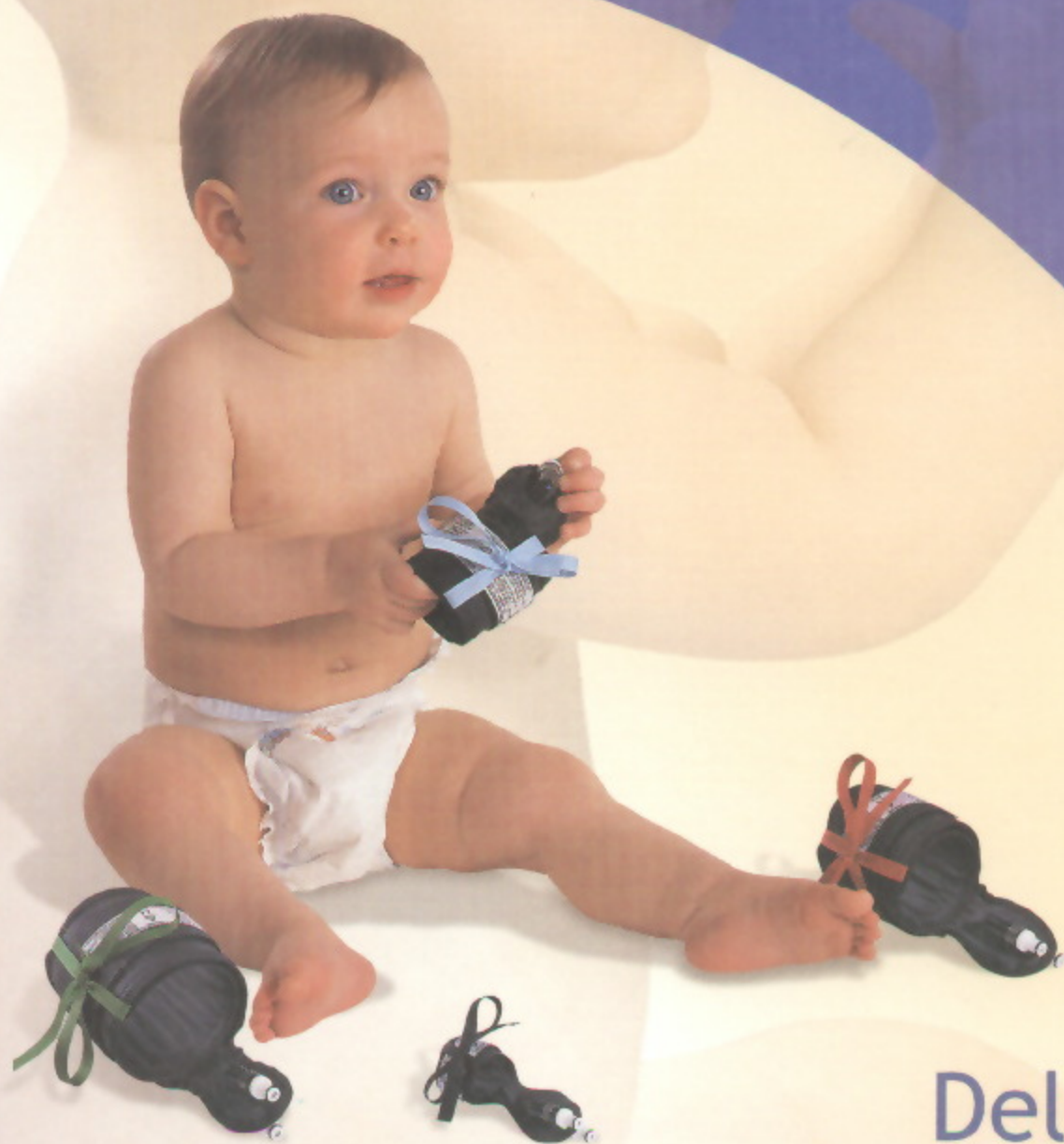

**Delfi** medical  
innovations inc.

# \* DELFI PEDIATRIC TOURNIQUET CUFFS WITH

CHILDREN ARE SPECIAL, AND THEY REQUIRE SPECIAL TOURNIQUETS FOR SURGERY.

DELFI'S PEDIATRIC TOURNIQUET CUFFS WITH MATCHING LIMB PROTECTION SLEEVES WERE INVENTED SPECIFICALLY FOR CHILDREN, WITH PATENTED FEATURES TO MEET THEIR UNIQUE REQUIREMENTS:

- FOUR SIZES OF PEDIATRIC TOURNIQUET CUFFS WITH MATCHING LIMB PROTECTION SLEEVES are designed to fit the limbs of most children, from neonates to adolescents
- LIMB PROTECTION SLEEVES THAT MATCH EACH CUFF SIZE are designed to help protect the skin and soft tissues from pinching and bruising
- P-SHAPED CUFF DESIGN reduces obstruction of pediatric surgical sites due to connectors and hoses
- WIDTH OF EACH CUFF is designed to reduce the need for high tourniquet pressures while minimizing obstruction of pediatric surgical sites
- SOFT AND THIN DESIGN of both Cuffs and Sleeves allows close conformance around pediatric limbs for the application of tourniquet pressure
- TWO INDEPENDENT CUFF FASTENERS make proper Cuff application easier and provide a secondary means of securing the Cuff around the limb during use
- PATENTED POSITIVE LOCKING CONNECTORS are fully compatible with A.T.S.\* Automatic Tourniquet Systems



Products covered by U.S. Patent No. 5,741,295; 5,649,954; 5,578,055; 5,454,831 and foreign patents and patent applications.

\* A.T.S. is a Registered Trademark of Zimmer Inc. Delfi Medical Innovations Inc. is not affiliated with Zimmer Inc.

# MATCHING LIMB PROTECTION SLEEVES \*

DELFI LIMB PROTECTION SLEEVES ARE DESIGNED TO HELP PROTECT A CHILD'S LIMB FROM POSSIBLE WRINKLING, PINCHING AND SHEARING OF SKIN AND SOFT TISSUES.

Each Delfi Limb Protection Sleeve is intended for use with a specific Delfi Pediatric Tourniquet Cuff.

Within the range of limb sizes recommended for the matching Cuff, the Sleeve helps to protect the limb's soft tissues beneath the Cuff while applying a low pressure that is less than a snugly wrapped Cuff.


The specially tapered end of the Sleeve is designed to help protect the limb from the Cuff tubing connectors.


After selecting the widest Delfi Pediatric Tourniquet Cuff that properly fits the limb and leaves a safe distance between the distal edge of the cuff and the joint, the correct Sleeve size is chosen by matching the color of the Sleeve to the color of the tie ribbon on the Cuff.

Each set of Delfi Pediatric Tourniquet Cuffs includes 4 Cuffs. The range of pediatric limb circumferences that each size of Cuff accommodates is shown in the Cuff Selection Chart.



CUFF SELECTION CHART

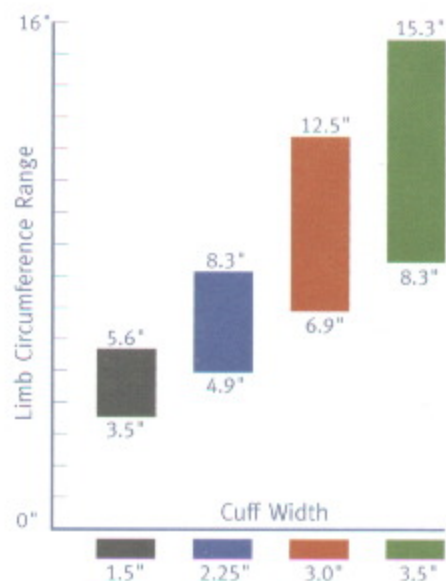

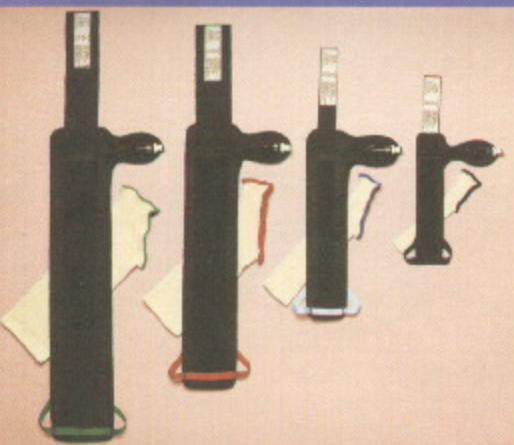

CATALOG NUMBER	DESCRIPTION
9-7800-005	Set of 4 DUAL PORT Pediatric Cuffs with Limb Protection Sleeves, Ports have Positive Locking Connectors
9-7700-005	Set of 4 SINGLE PORT Pediatric Cuffs with Limb Protection Sleeves, Ports have Positive Locking Connectors
9-7900-001	Limb Protection Sleeve P-1.5 (Black) Box of 10
9-7900-002	Limb Protection Sleeve P-2.25 (Blue) Box of 10
9-7900-003	Limb Protection Sleeve P-3.0 (Red) Box of 10
9-7900-004	Limb Protection Sleeve P-3.5 (Green) Box of 10

Phone: 1-800-933-3022 (Toll Free in U.S. and Canada)  
 +604-742-0600 (Elsewhere)  
 Fax: +604-742-3800 (Global)  
 Email: info@delfimedical.com  
 Internet: www.delfimedical.com  
 Address: 106-1099 W. 8th Ave. Vancouver, Canada V6H 1C3

**Delfi** medical  
 innovations inc.



DESIGNED IN CANADA M5375